

Data Sheet

ATAD2A Inhibitor Screening Assay Kit

Catalog # 32601

DESCRIPTION: The *ATAD2A Inhibitor Screening Assay Kit* is designed to measure the inhibition of ATAD2A binding to its substrate. The *ATAD2A Inhibitor Screening Assay Kit* comes in a convenient AlphaLISA[®] format, with biotinylated histone peptide substrate, assay buffer, detection buffer and purified GST-tagged ATAD2A bromodomain to perform a total of 384 enzyme reactions. The key to the *ATAD2A Inhibitor Screening Assay Kit* is a specific binding of the ATAD2A bromodomain to the acetylated histone substrate. With this kit, only three simple steps on a microtiter plate are required. First, a sample containing ATAD2A and an inhibitor of choice is incubated with the biotinylated substrate for thirty minutes. Next, acceptor beads are added, then donor beads, followed by reading the Alpha-counts.

COMPONENTS:

Catalog #	Component	Amount	Storage	
31121	GST-ATAD2A (981-1108)	200 µg	-80 °C	(Avoid freeze/thaw cycles!)
	BET Bromodomain Ligand	400 µl	-80 °C	
	Non-acetylated Ligand 1	200 µl	-80 °C	
	3x ATAD2A assay buffer	4 ml	-20 °C	
	3x ATAD2A detection buffer	3 ml	-20 °C	

MATERIALS OR INSTRUMENTS REQUIRED BUT NOT SUPPLIED:

Glutathione AlphaLISA[®] Acceptor Beads, 5 mg/ml (PerkinElmer #AL109C)
AlphaScreen[®] Streptavidin-conjugated donor beads, 5 mg/ml (PerkinElmer #6760002S)
Optiplate -384 (PerkinElmer #6007290)
AlphaScreen[®] microplate reader
Adjustable micropipettor and sterile tips

APPLICATIONS: Useful for the study of bromodomain binding assays, screening inhibitors and selectivity profiling.

CONTRAINDICATIONS: Green and blue dyes that absorb light in the AlphaScreen[®] signal emission range (520-620 nm), such as Trypan Blue. Avoid the use of the potent singlet oxygen quenchers such as sodium azide (NaN₃) or metal ions (Fe²⁺, Fe³⁺, Cu²⁺, Zn²⁺ and Ni²⁺). The presence of >1% RPMI 1640 culture medium leads to a signal reduction due to the presence of excess biotin and iron in this medium. MEM, which lacks these components, does not affect AlphaScreen[®] assays.

STABILITY: At least one year from date of receipt when stored as directed.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082**, Fax **1.858.481.8694**
Or you can Email us at: info@bpsbioscience.com
Please visit our website at: www.bpsbioscience.com

REFERENCE: Filippakopoulos, P., *et al.*, *Cell* 2012; **149**:214.

ASSAY PROTOCOL:

All samples and controls should be tested in duplicate. Use slow shaking for all incubations.

Step 1:

- 1) Prepare the master mixture: N wells × (2.5 µl **3x ATAD2A assay buffer** + 1 µl **BET Bromodomain Ligand** + 1.5 µl **H₂O**).
- 2) Thaw **ATAD2A** on ice. Upon first thaw, briefly spin tube containing protein to recover full content of the tube. Aliquot protein into single use aliquots. Store remaining undiluted protein in aliquots at -80°C immediately. *Note: ATAD2A is very sensitive to freeze/thaw cycles. Do not re-use thawed aliquots or diluted protein.*
- 3) Dilute **ATAD2A** in 1x **ATAD2A assay buffer** at 200 ng/µl. Keep diluted protein on ice until use. Discard any unused diluted protein after use.
- 4) Add 5 µl of master mixture to each well designated for the “Positive Control”, “Test Inhibitor”, and “Blank”. For the “Substrate Control”, add 2.5 µl **3x ATAD2A assay buffer** + 1.5 µl **H₂O** + 1 µl of **Non-acetylated Ligand 1**.

	Blank	Substrate Control	Positive Control	Test Inhibitor
3x ATAD2A assay buffer	2.5 µl	2.5 µl	2.5 µl	2.5 µl
BET Bromodomain Ligand	1 µl	–	1 µl	1 µl
Non-acetylated Ligand 1	–	1 µl	–	–
H ₂ O	1.5 µl	1.5 µl	1.5 µl	1.5 µl
Test Inhibitor/Activator	–	–	–	2.5 µl
Inhibitor buffer (no inhibitor)	2.5 µl	2.5 µl	2.5 µl	–
1x ATAD2A assay buffer	2.5 µl			
ATAD2A (200 ng/µl)	–	2.5 µl	2.5 µl	2.5 µl
Total	10 µl	10 µl	10 µl	10 µl

- 5) Add 2.5 µl of **inhibitor solution** to each well designated “Test Inhibitor”. For the “Positive Control”, “Substrate Control” and “Blank”, add 2.5 µl of the same **solution without inhibitor** (inhibitor buffer). *Note: Keep DMSO concentration below 0.5 %.*
- 6) Add 2.5 µl of **1x ATAD2A assay buffer** to the well designated “Blank”.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082**, Fax **1.858.481.8694**
 Or you can Email us at: info@bpsbioscience.com
 Please visit our website at: www.bpsbioscience.com

- 7) Initiate reaction by adding 2.5 μ l of diluted **ATAD2A** prepared as described above. Incubate at room temperature for 30 minutes.

Step 2:

Note: Protect your samples from direct exposure to light!


- 1) Dilute Glutathione AlphaLISA Acceptor Beads (PerkinElmer #AL109C) 250-fold with 1x ATAD2A detection buffer. Add 10 μ l per well. Shake plate briefly. Incubate at room temperature for 30 minutes.

Step 3:

- 1) Dilute Streptavidin-conjugated donor beads (PE #6760002S) 125-fold with 1x ATAD2A detection buffer. Add 10 μ l per well. Incubate at room temperature for 15 - 30 minutes.
2) Read Alpha-counts.

Due to lot to lot variability in AlphaScreen[®] bead performance, it may be necessary to optimize assay conditions. For example, slight adjustments to bromodomain or ligand concentrations may improve signal-to-noise ratio.

EXAMPLE OF ASSAY RESULTS:

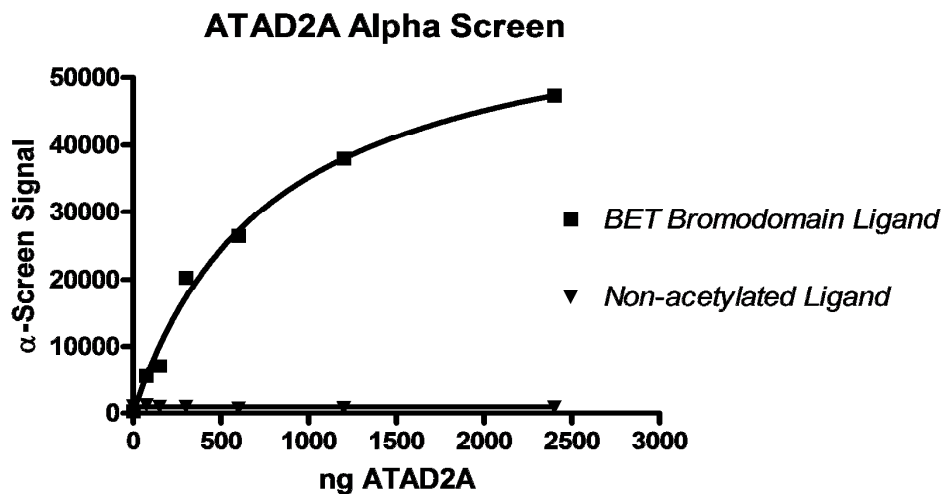

ATAD2A binding activity, measured using the ATAD2A Inhibitor Screening Assay Kit, BPS Bioscience, Catalog # 32601. *Data shown is lot-specific. For lot-specific information, please contact BPS Bioscience, Inc. at info@bpsbioscience.com.*

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082**, Fax **1.858.481.8694**

Or you can Email us at: info@bpsbioscience.com

Please visit our website at: www.bpsbioscience.com



6044 Cornerstone Ct. West, Ste. E
San Diego, CA 92121
Tel: 1.858.829.3082
Fax: 1.858.481.8694
Email: info@bpsbioscience.com

RELATED PRODUCTS:

<u>Product Name</u>	<u>Catalog #</u>	<u>Size</u>
ATAD2A (981 – 1108), His-tag	31109	100 µg
ATAD2A (981 – 1108), GST-tag	31121	100 µg
ATAD2B (953 – 1080), His-tag	31117	100 µg
BAZ2B (2054 – 2168), His-tag	31113	100 µg
BRD2 (339 – 459), His-tag*	31020	100 µg
BRD3 (29 – 145), His-tag*	31030	100 µg
BRD3 (306 – 417), His-tag*	31031	100 µg
BRD4 (49 – 170), His-tag*	31042	100 µg
BRD4 (342 – 460), His-tag*	31043	100 µg
BRD9 (135 – 242), His-tag	31090	100 µg
BRDT (22 – 138), His-tag*	31101	100 µg
BRDT (257 – 382), His-tag	31100	100 µg
BET Bromodomain Ligand	33000	0.5 mL
BRD2 (BD2) Inhibitor Screening Kit	32522	384 rxns.
BRD4 (BD1) Inhibitor Screening Kit	32514	384 rxns.
BRD4 (BD2) Inhibitor Screening Kit	32524	384 rxns.
BRD3 (BD1) Inhibitor Screening Kit	32513	384 rxns.
BRD3 (BD2) Inhibitor Screening Kit	32523	384 rxns.
BAZ2B Inhibitor Screening Kit	32600	384 rxns.
(+)-JQ1 Bromodomain Inhibitor	27400	10 mg

*Also available with GST-tag

AlphaScreen® and AlphaLISA® are registered trademarks of PerkinElmer, Inc.

OUR PRODUCTS ARE FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

To place your order, please contact us by Phone **1.858.829.3082**, Fax **1.858.481.8694**

Or you can Email us at: info@bpsbioscience.com

Please visit our website at: www.bpsbioscience.com