



To the bottom of tubes!

Super Long Tip

For more information : http://www.funakoshi.co.jp/exports_contents/80530

Super Long Tip is a long pipette tip with 15.5 cm. It is attachable to most suppliers' micropipettes (100 µL - 300 µL size) and repeating pipettes, such as eppendorf's combitips plus®.

Features

- Long tip - can reach the bottom of centrifuge tube (15 mL and 50 mL). 15.5 cm = 6.1 inch.
- Compatible with most micropipettes (100 µL - 300 µL size), repeat pipettes, such as eppendorf's combitips plus®.
- Available in two formats - rack and bag type.
- Both tip and rack are autoclavable.



Real Scale



Product Information

[Manufacturer : NIP]

Product Name	Size	Code
Super Long Tip 15.5cm, Sterilized	200 pcs/bag	891141
	96 pcs/rack	891143
Super Long Tip 15.5cm, Non-Sterile	200 pcs/bag	891140
	96 pcs/rack	891142

Ice On

For more information : http://www.funakoshi.co.jp/exports_contents/80502

Ice On is a sample tube cooling rack, made of Aluminium. As tube walls adhere to rack tightly, tubes are cooled more effectively than putting tubes onto ice directly.


Features

- Higher cooling efficiency than regular on ice cooling.
- Surface-treated for avoiding corrosion.
- Compact and Stackable




Principle

Regular cooling on ice

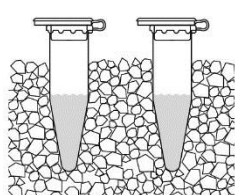

Gap between ice and tube makes cooling efficiency low.

Cooling by using Ice On

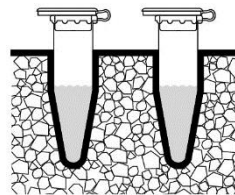

Higher cooling efficiency
 • Great thermal conductivity of cold by aluminium.
 • No gap between tube and rack – rapid cooling

Product	Ice On Type I (# IO-1)	Ice On Type II (# IO-2)
Material	Aluminum Die-cast	
Storage	1.5 ml x 24	0.2 ml x 96 (0.2 ml 8-strip x 12, 0.2 ml 12-strip x 8 or 0.2 ml 96well plate x 1) , 0.5 ml x 6 or 1.5 ml x 4
Size	169 x 69 x 47 mm	177 x 95 x 47 mm
Weight	300 g	470 g

Product Information

[Manufacturer : SKI]

Product Name	Size	Code
Ice On Type I	1 piece	IO-1
Ice On Type II	1 piece	IO-2

NOTE

* All products on here are for research use only, not for diagnostic use.
 * Products might be changed for improvement without notice.

* All brand names and product names are registered trademarks of their respective companies.

* Please contact your local distributors for orders, quote request and any inquiry.

Your Local Distributor

Funakoshi Co., Ltd.

Address: 9-7 Hongo 2-Chome, Bunkyo-ku,
 Tokyo 113-0033 JAPAN
 Phone : +81-3-5684-6296
 Fax : +81-3-5684-6297
 Email : export@funakoshi.co.jp

FNA-1801-F03